

**Fig 2.** Schematic drawing of MicroCal iTC<sub>200</sub>.

With MicroCal Auto-iTC<sub>200</sub>, all filling, injection, and cell-cleaning functions are fully automated, controlled, and operated through software that includes integrated experiment design wizards to assist in selecting experimental parameters.

Data analysis is performed with Origin® software using fitting models to calculate reaction stoichiometry ( $n$ ), binding constant ( $K_D$ ), enthalpy ( $\Delta H$ ), and entropy ( $\Delta S$ ), and the results are presented in Microsoft™ Excel™ format for further analysis or data transfer.


## The ITC experiment

ITC measures the heat absorbed or generated when molecules interact. These heat changes are small, typically submillionths of a degree, but are universal and can be detected by these very sensitive instruments.


A solution of the biomolecule is first placed in the sample cell and a ligand solution in a matching buffer is placed in the syringe (Fig 3A). When the ligand solution is injected into the cell, the ITC instrument detects heat that is released or absorbed as a result of the interaction. This is done by measuring the changes in the power needed to maintain isothermal conditions between the reference and the sample cell. Injections are performed repeatedly, and result in peaks that become smaller as the biomolecule becomes saturated. Eventually, the peak sizes remain constant and represent only the heats of dilution.

Once titration is completed, the individual peaks are integrated by the instrument software (Fig 3B) and presented in a Wiseman plot (Fig 3C). An appropriate binding model is chosen and the isotherm is fitted to yield the binding enthalpy  $\Delta H$ , the  $K_D$ , and the stoichiometry,  $n$ . From these data, Gibbs' free energy,  $\Delta G$  and entropy,  $\Delta S$  are calculated.

Besides confirming direct interactions with the target of interest, thermodynamic measurements also provide insight into the nature of the noncovalent forces responsible for binding. Polar interactions tend to contribute favorably to the enthalpic component, whereas entropically favored interactions tend to be more hydrophobic. Figure 4 shows representative ITC binding isotherms for two interactions with the same affinity but with different mechanisms of binding.

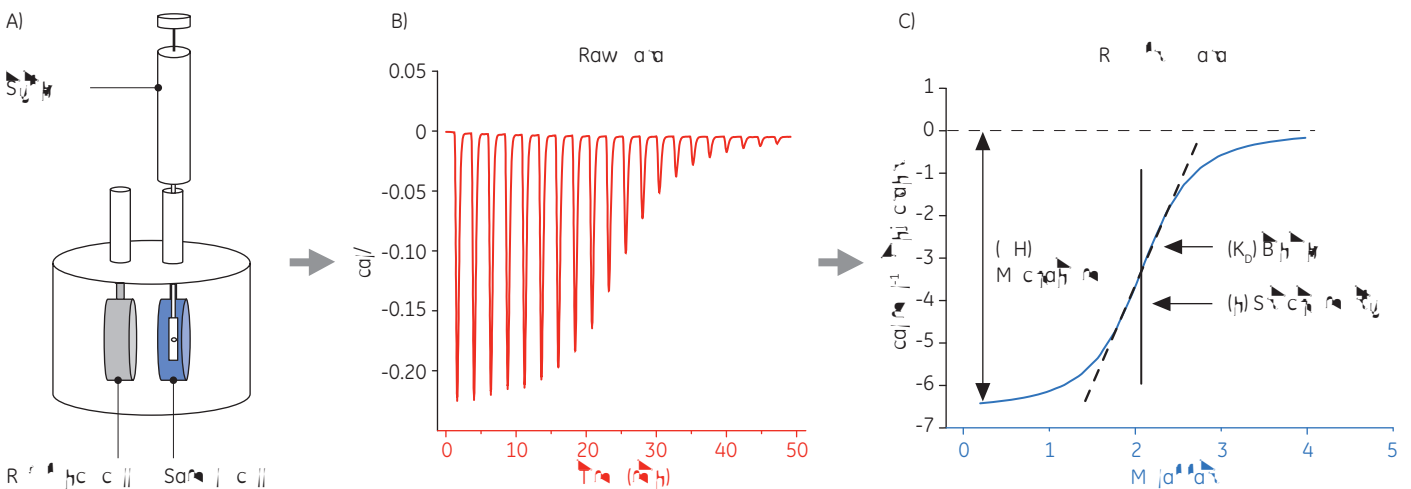

**Fig 3.** (A) The ligand is titrated into the sample cell. (B) An exothermic reaction releases heat and gives negative peaks. (C) The peaks are integrated and presented in a Wiseman plot.

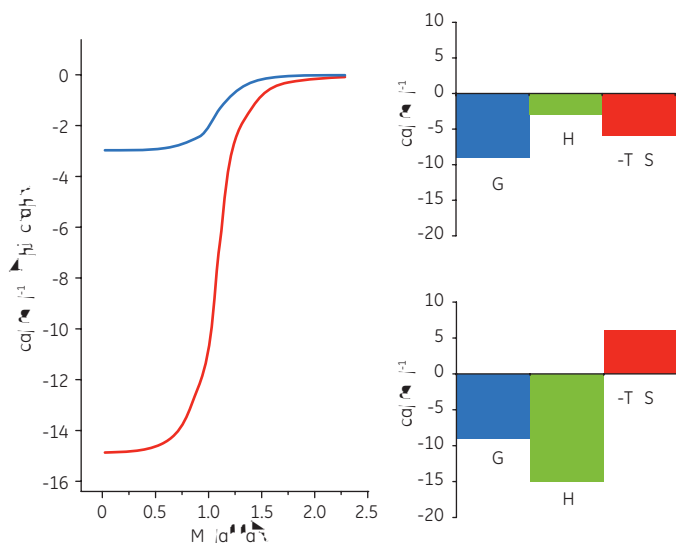


**Fig 4.** An entropy-driven interaction (blue curve and top right profile) tends to be more hydrophobic in character compared to an enthalpy-driven interaction (red curve and bottom right profile), which tends to be driven by hydrogen-bonding and van der Waals interactions.

## Elucidation of reaction mechanisms


To demonstrate how ITC can be used not only to measure binding affinity but also provide insight into the mechanism of binding, the interaction of peptides with the protein target Bcl-2 (name derived from B-cell lymphoma 2) has been studied with MicroCal iTC<sub>200</sub>.

Experiments were carried out at 25°C. The sample cell was filled with Bcl-2 (30 μM solution) in 50 mM HEPES, pH 7.4, 100 mM NaCl, 0.5mM TCEP, and 5% DMSO. The peptides were diluted to a concentration of 250 μM in the same buffer. The injection volumes were 3 μl each, injection time 6 s, and a 150 s delay between each injection. The results are shown in Figure 5.

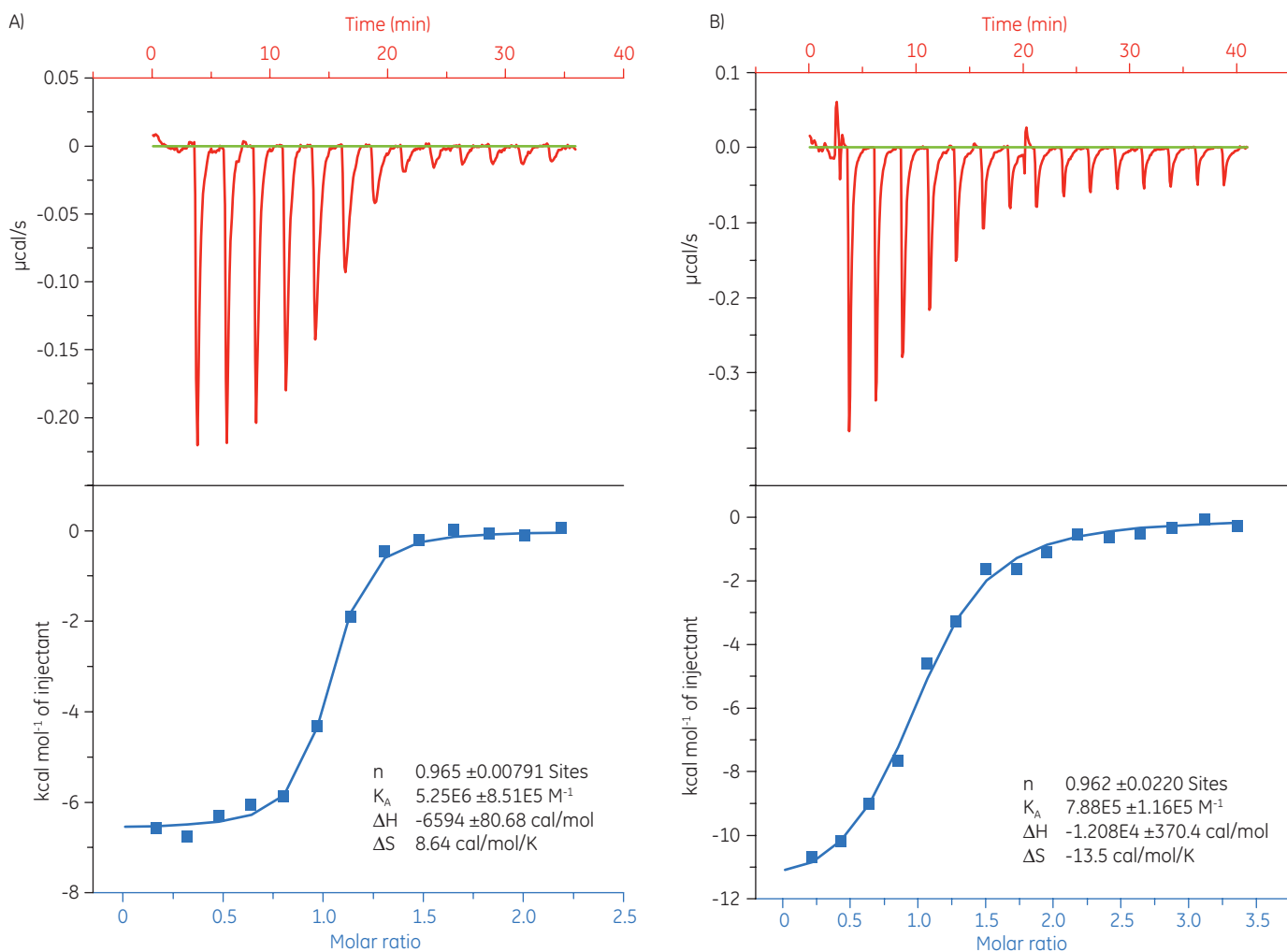

**Fig 5.** The top panels show the raw data. The bottom panels show the binding isotherms created by plotting the integrated heat peaks against the molar ratio of the peptide. (A) Injection of BAD-like peptide. (B) Injection of BAX peptide.

The binding affinity of the BAD-like (Bcl-2-associated death promoter) peptide for Bcl-2 protein is approximately six-fold stronger than that of the BAX (Bcl-2-associated X protein) peptide. Visualization of thermodynamic parameters in the form of a binding signature plot (Fig 6) makes it easier to see how the enthalpic and entropic components contribute to the overall affinity, represented here by  $\Delta G$ . These plots reveal that the binding of the BAD-like peptide to Bcl-2 is comprised of hydrogen bonding and hydrophobic interactions as indicated by the negative or favorable binding enthalpy ( $\Delta H$ ) and entropy factor ( $T\Delta S$ ), while the binding of BAX involves more conformational changes as indicated by the unfavorable entropy in addition to hydrogen bonding ( $\Delta H$ ).

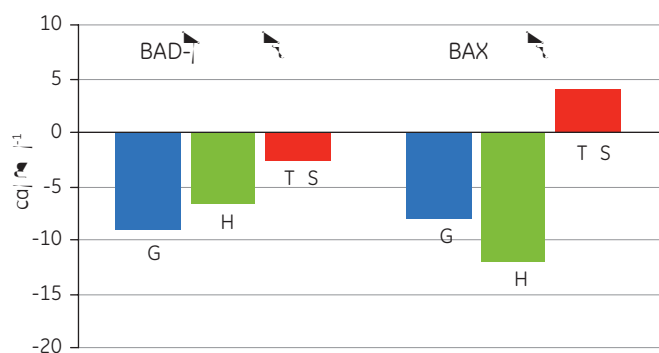

**Fig 6.** The binding signature (free energy, binding enthalpy, and entropy factor) plotted for the two binding events.


## Simultaneous determination of thermodynamic parameters

To study the interaction of carbonic anhydrase with five known inhibitors, a series of 20 runs was performed using a MicroCal Auto-iTC<sub>200</sub> (Table 1). The workflow, including sample introduction, titration, and cleaning was completely automated. Each inhibitor was run four times and all 20 runs were completed in 15 h. The excellent reproducibility of the replicate injections is demonstrated in the overlaid binding isotherms (Fig 7).


Both CBS and furosemide have similar binding affinities, yet different enthalpies, suggesting that the binding mechanisms are different.

**Table 1.** Thermodynamic parameters determined for the interaction of five inhibitors with bovine carbonic anhydrase II (BCA). The values are the average of four separate runs with errors shown

Ligand/ Concentration	BCA ( $\mu\text{M}$ )	n	$K_D$ ( $\mu\text{M}$ )	$\Delta G$ ( $\text{kcal mol}^{-1}$ )	$\Delta H$ ( $\text{kcal mol}^{-1}$ )	$-T\Delta S$ ( $\text{kcal mol}^{-1}$ )
Acetazolamide/ 0.126 mM	10	$0.98 \pm 0.02$	0.06	-9.87	$-11.15 \pm 0.46$	1.28
CBS/ 0.414 mM	30	$1.00 \pm 0.04$	0.96	-8.21	$-10.19 \pm 0.12$	1.98
Furosemide/ 0.426 mM	30	$0.98 \pm 0.08$	0.92	-8.23	$-7.06 \pm 0.20$	-1.17
Sulfanilimide/ 0.441 mM	30	$0.99 \pm 0.05$	4	-7.35	$-7.93 \pm 0.39$	0.58
TFMSA/ 0.525 mM	30	$1.03 \pm 0.02$	0.35	-8.8	$-2.03 \pm 0.07$	-6.77

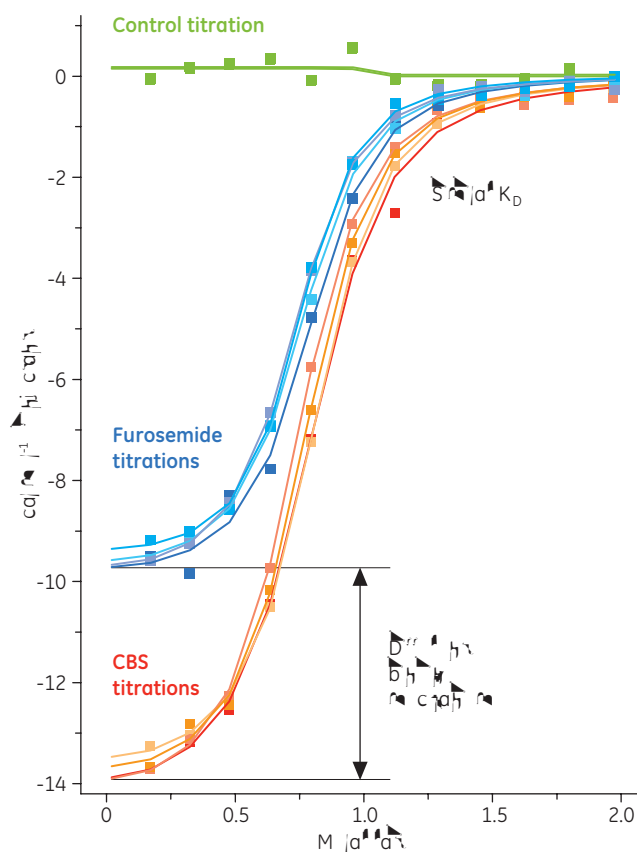

**Fig 7.** Overlay plots from four repeated titrations of bovine carbonic anhydrase II, with CBS (left) and furosemide (right).

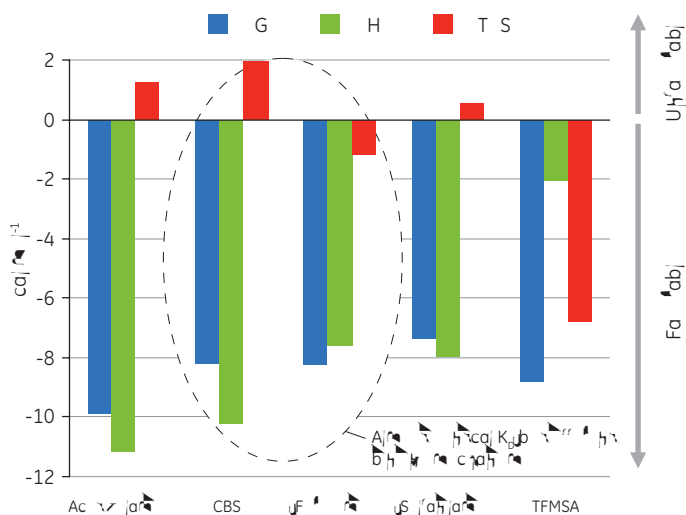


**Fig 8.** Binding signatures for five inhibitors of bovine carbonic anhydrase II.

The binding signatures for these interactions are shown in Figure 8. From the enthalpy and unfavorable entropy factor, CBS has binding based entirely on hydrogen and van der Waals bonds with some conformational changes reducing the affinity. Furosemide has a more balanced binding based on hydrogen and van der Waals bonds as well as hydrophobic effects. This illustrates that even when the inhibitors have almost identical binding affinities, important differences in binding mechanisms can be easily determined.


## Verify hits through increased understanding of binding mechanisms

It is important to rule out false positives from a primary screening at an early stage. ITC can provide such information, saving considerable time and effort later on. A 20  $\mu\text{M}$  solution of target protein (TP) was loaded into the sample cell of MicroCal iTC<sub>200</sub> and titrated with Compound X (Fig 9). The  $K_D$  was determined to 4.9  $\mu\text{M}$ , which correlated well with studies performed with Biacore™ systems and NMR, thus confirming that Compound X is suitable for further studies.

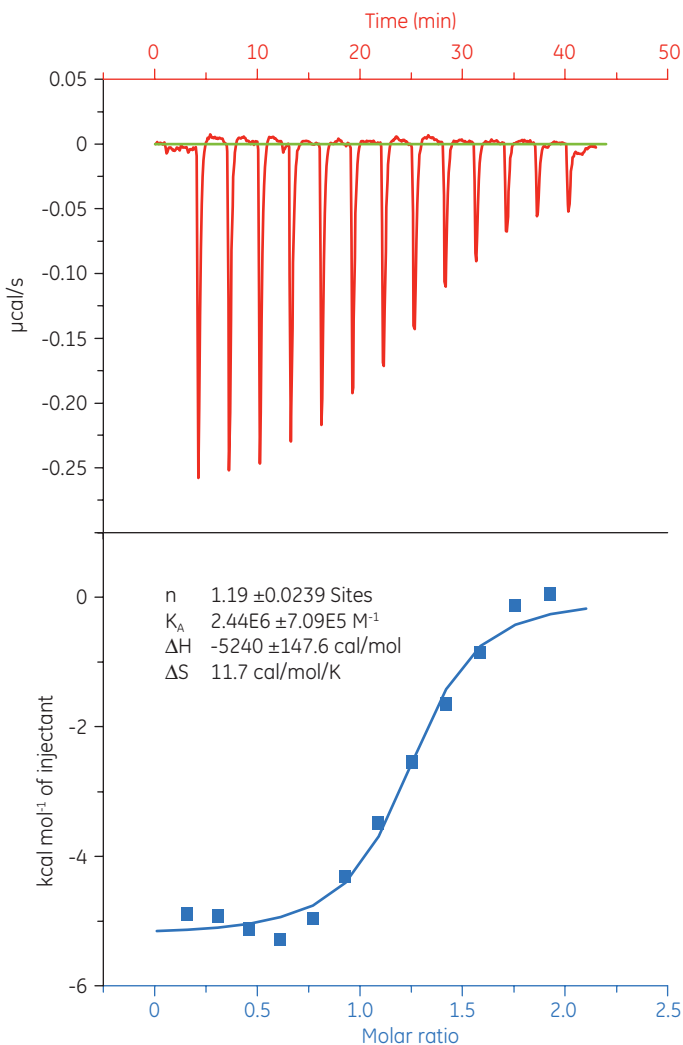

**Fig 9.** Raw data and binding isotherm for the interaction of Compound X with TP, data show a reversible binding.


When the same target solution was titrated with Compound Y, the results were very different (Fig 10). In the left panel, Compound Y was titrated with TP. The isotherm shows a binding affinity of 120 nM, but the binding enthalpy was about 1000-fold larger than expected and the stoichiometry value (0.01) very low. In the right panel, the same drug candidate was titrated with bovine serum albumin (BSA). Taken together, the results indicate nonspecific activity and this compound is not suitable for further development. Based on these experiments, Compound Y was considered to be a false positive and was rejected for further study.

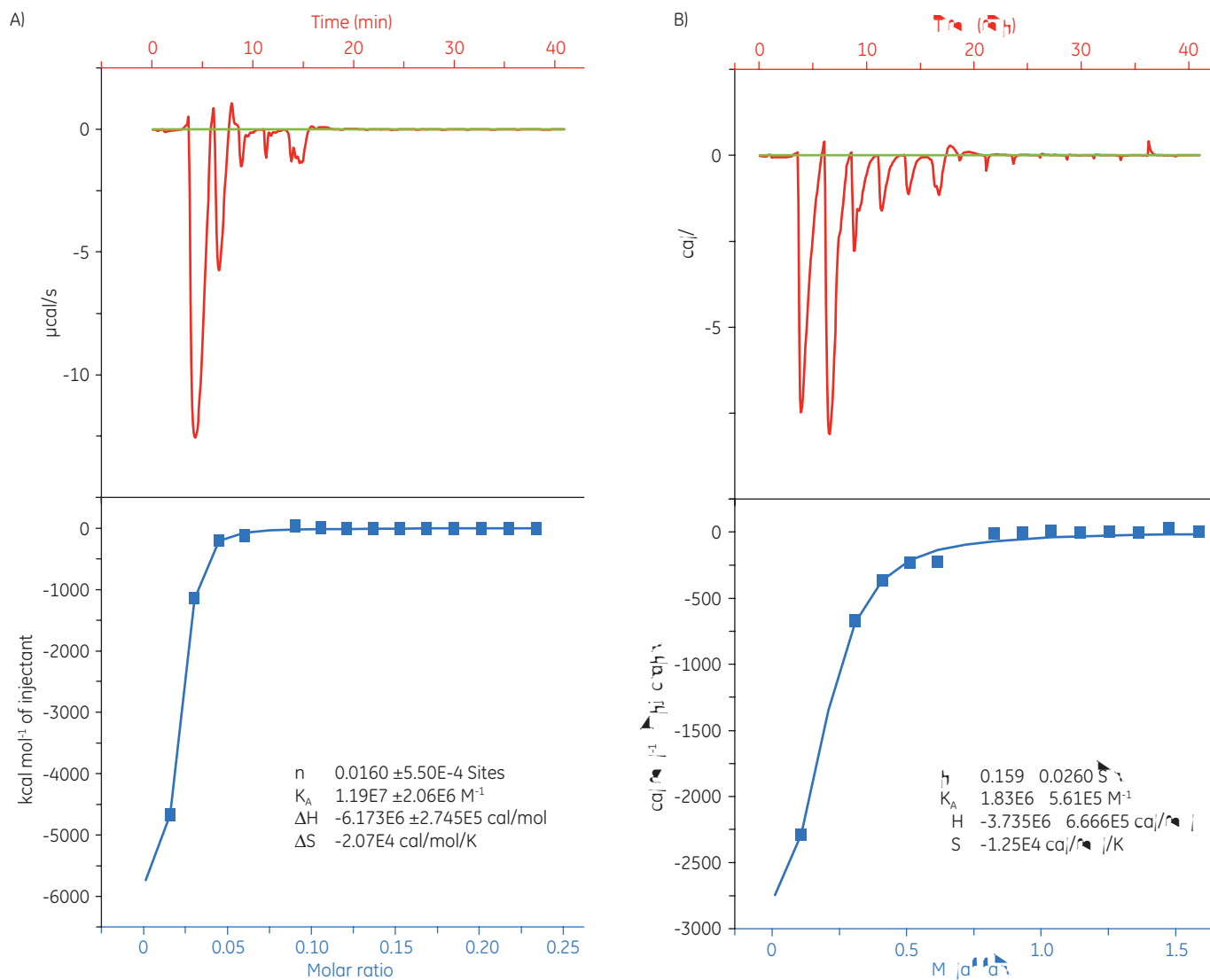

**Fig 10.** Data for an irreversible interaction of Compound Y with TP (left) and nonspecific binding to BSA (right).

## Assessing protein quality

ITC can be used to assess the activity level of a target protein before use in high throughput screening. Two batches of a target protein were compared using MicroCal iTC<sub>200</sub> by titration with a standard peptide (Fig 11).

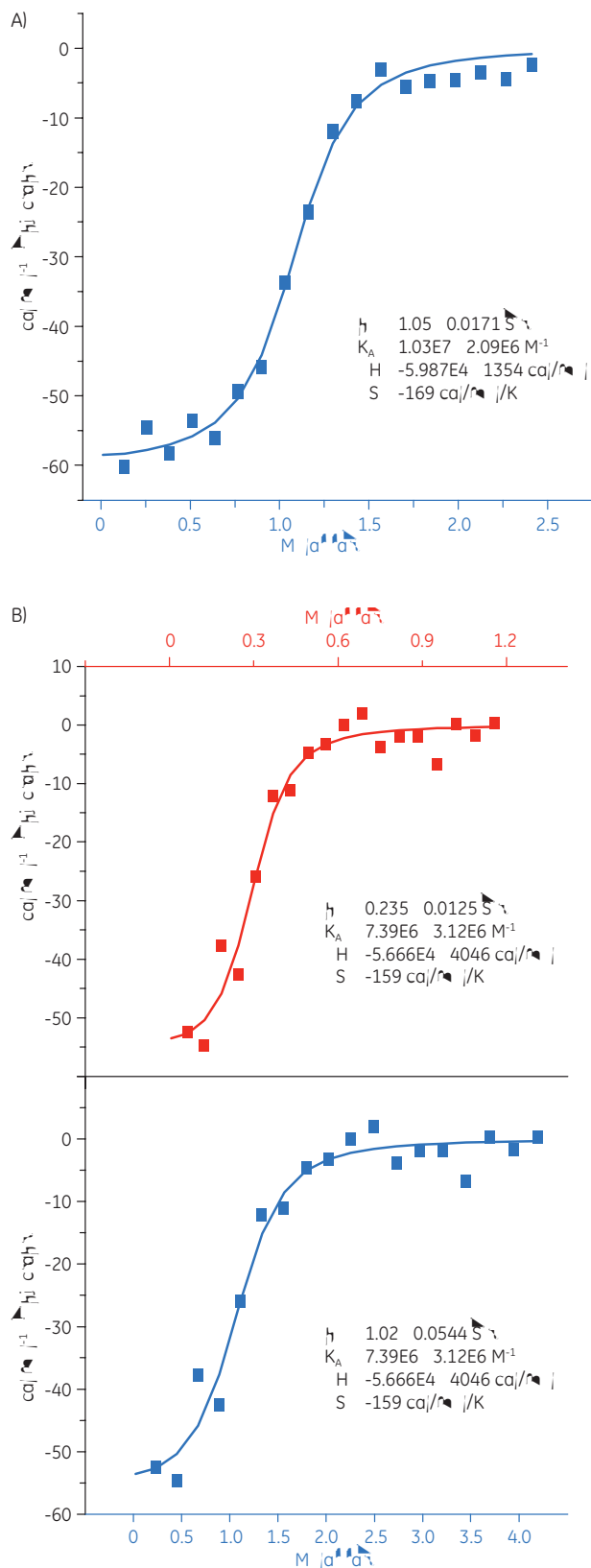

**Fig 11.** Two batches of target protein (A and B, respectively) titrated with a standard peptide. The sample cell contained the protein at a 10  $\mu$ M concentration, and the peptide solution was 50  $\mu$ M.

The result for Batch 1 (Fig 11 A) demonstrates a typical isotherm with a  $K_D$  of 97 nM and  $n = 1$ . The second batch (Fig 11 B) has a  $K_D$  of 135 nM but  $n$  is only 0.23. The analysis of the same set of data but with an estimated protein concentration of 2.3  $\mu$ M gives the same  $K_D$  value and  $n = 1$ . This indicates that 75% of the Batch 2 protein was inactive. The Batch 2 protein was rejected for use in the screening.

## Acknowledgement

The data shown in Figs 5, 6, 9, and 10 were kindly provided by Dr. Lin Gao, Hoffman-La Roche, Nutley NJ, USA.

## Technical specifications

### MicroCal Auto-iTC<sub>200</sub> and MicroCal iTC<sub>200</sub>

Sample volume	
MicroCal Auto-iTC <sub>200</sub>	370 $\mu$ L
MicroCal iTC <sub>200</sub>	280 $\mu$ L
Equilibration time from 25° to 5°C	< 6 min
Response time	10 s*
Throughput	
MicroCal Auto- iTC <sub>200</sub>	Up to 75 per 24 h (SIM)
MicroCal iTC <sub>200</sub>	8-12 per 8 h day
Injection syringe volume	40 $\mu$ L
Cell material	Hastelloy
Cell configuration	Coin-shaped
Cell volume	200 $\mu$ L
Noise	0.2 ncal/s†
Temperature range	2°C to 80°C
Temperature stability	150 $\mu$ cal at 25°C
Multiple feedback modes	Yes (passive, high gain, low gain)
Automated upgrade available	
MicroCal Auto- iTC <sub>200</sub>	NA
MicroCal iTC <sub>200</sub>	Yes
Weight	
MicroCal iTC <sub>200</sub>	9.4 kg
MicroCal Auto-iTC <sub>200</sub>	90.7 kg
Dimensions (W × H × D)	
MicroCal iTC <sub>200</sub>	21.0 × 33.7 × 34.9 cm
MicroCal Auto-iTC <sub>200</sub>	62.5 × 76.8 × 57.2 cm

\* high feedback

† cell temperature = 25°C, no feedback, 5 s filter period, stir speed = 1000 rpm

## Ordering information

Product	Code number
MicroCal iTC <sub>200</sub> system, controller and software	28-4289-55
MicroCal Auto-iTC <sub>200</sub> system, controller and software	28-4289-83
MicroCal iTC <sub>200</sub> upgrade to MicroCal Auto-iTC <sub>200</sub>	28-4289-84